THE LED 2.0 BIOHACK MASK REPRESENTS A NEW APPROACH TO SKINCARE BASED ON THE PRINCIPLES OF LONGEVITY, AS IT SUPPORTS REGENERATION, VITALITY, AND THE LONG-TERM QUALITY OF THE SKIN.

Envy Therapy

LED 2.0 Biohack Mask

RADICAL
SKINCARE SYSTEM

[EN]

This instruction manual contains an explanation of the individual functions of the mask and warnings regarding its use.

Before using the mask, carefully read this instruction manual and familiarize yourself with its contents.

Charge the mask before first use. The mask is not a medical device, but a cosmetic tool.

The effect of the mask is influenced by skin type and the type of products you use.

1. PRODUCT INFORMATION
 - 1.1. MAIN BENEFITS
 - 1.2. MODE DESCRIPTION
2. USAGE
 - 2.1. HOW TO USE THE MASK
 - 2.2. MASK CONTROLS AND REMOTE CONTROL
 - 2.3. EXPOSURE TIME AND DURATION OF ONE THERAPY SESSION
 - 2.4. FREQUENCY OF USE
3. CONTRAINDICATIONS
4. WARNINGS FOR USE
5. PACKAGE CONTENTS
6. PRODUCT PARAMETERS

1. PRODUCT INFORMATION

The LED 2.0 Biohack Mask represents a technological shift in home skincare and a new approach based on the principles of longevity. The mask features an advanced infrared spectrum of 1070 nm, which penetrates significantly deeper into the skin than classic LED therapy and enables more intensive stimulation of regenerative processes in the deeper layers of the skin.

The mask combines the infrared spectrum with other precisely targeted light modes. It offers a comprehensive solution for signs of aging, pigmentation, and acne.

The mask is a cosmetic tool intended to support skin quality, not a medical device. Its effect depends on skin type, regularity of use, and a proper skincare routine.

1.1. MAIN BENEFITS

- Deeper light penetration (up to 1070 nm).
- 3 targeted modes according to the skin's needs.
- Wireless use.
- Comfortable fit with one adjustable strap.
- No mouth cutout for greater comfort.
- Freedom of movement during therapy.
- Maximum comfort and easy use.

1.2. MODE DESCRIPTION

The LED 2.0 Biohack Mask uses a 1070 nm infrared spectrum in combination with three targeted light modes:

Mode 1: Anti-Age (1070 nm + red light 630 nm)

Anti-Age mode for firmer and younger-looking skin. Stimulates collagen production, supports elasticity, and helps reduce fine lines. The skin gradually becomes more compact, more supple, and visibly smoother.


Mode 2: Anti-Acne (1070 nm + blue light 465 nm)

Anti-Acne mode regulates sebum production, acts against acne-causing bacteria, and helps alleviate inflammatory symptoms. An ideal solution for problematic or oily skin prone to imperfections.

Mode 3: Anti-Pigment (1070 nm + yellow light 590 nm)

Anti-Pigment mode is focused on pigmentation and uneven skin tone. It helps reduce redness, promotes brightening, and evens out the skin. Ideal for tired, dull, or pigmented skin that needs to balance its tone and regain freshness.

Infrared spectrum technology with penetration of up to 1070 nm allows deeper light penetration into the skin and activation of regenerative processes in the deeper layers of the skin, thereby supporting the principles of longevity skincare. The combination of precisely targeted modes and the synergistic action with active ingredients from a booster or serum creates a new standard of home skin biohacking. LED light is safe, non-invasive, and suitable for home use. Thanks to advanced technology, it can also work in deeper layers of the skin, which are typically the target of laser technologies.

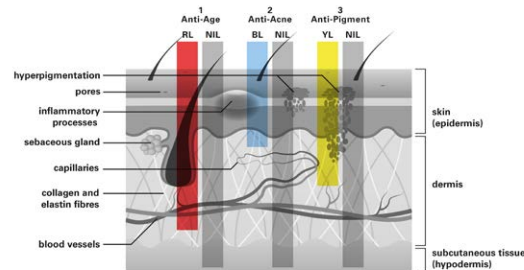BL (Blue light) 465 nm
YL (Yellow light) 590 nm

RL (Red light) 630 nm
NIL (NIR-Infrared light) 1070 nm

2. USAGE

2.1. HOW TO USE THE MASK

1. Thoroughly remove makeup and cleanse the skin.

We recommend using a gentle cleansing product Envy Therapy Cleansing Face Milk, which removes makeup, sebum, and impurities without disrupting the skin barrier.

2. Apply a serum or booster.

Before putting on the mask, apply a lightweight, non-greasy serum or concentrated booster suitable for light therapy. Active ingredients applied before therapy work more intensively because the infrared spectrum supports microcirculation and cellular metabolism, thereby increasing their effectiveness.

Suitable active ingredients: Hyaluronic acid, niacinamide, peptides, ceramides, panthenol (vitamin B5), aloe vera, beta-glucan, centella asiatica, allantoin, growth factors, exosomes, and NAD⁺, especially hydrating, regenerating, and soothing active ingredients.

Unsuitable active ingredients: Retinoids (retinol, tretinoin, adapalene, tazarotene), AHA/BHA/PHA acids (e.g. glycolic acid, lactic acid, mandelic acid, salicylic acid, gluconolactone, lactobionic acid), azelaic acid, strong forms of vitamin C (ascorbic acid), hydroquinone, and all chemical, mechanical, as well as enzymatic peels.

3. Put on the mask and secure it using the adjustable strap.

The mask must fit comfortably against the face so that the light is distributed evenly.

4. Turn on the mask using the remote control or the controls directly on the mask and select the desired mode.

Choose the mode according to the skin's current needs.

5. During therapy, you may move freely and let the active ingredients work.

Light energy supports cellular regeneration and at the same time enhances the effect of the applied product. The combination of light and active ingredients ensures a more intensive and targeted result.

6. After completing the therapy, apply a moisturizing cream.

Lock in the effect of the therapy with a suitable cream that supports hydration and restoration of the skin barrier.

2.2. MASK CONTROLS AND REMOTE CONTROL

The mask can be controlled in two ways: by using the remote control or directly via the button built into the mask.

Controls directly on the mask

There is an integrated control button located on the right side of the mask.

- Turning on the mask: press and hold the button for approximately 2–3 seconds.
- Switching modes: briefly pressing the button switches between the individual light modes (Mode 1, Mode 2, Mode 3).
- Turning off the mask: press and hold the button again for a few seconds until the device switches off.

This system enables simple operation even without using the remote control.


Remote control

The mask is equipped with a separate remote control. The remote control allows switching between the three light modes (MODE button), but cannot turn the mask on and off. You must always turn it on and off using the integrated button on the mask.

MODE button
Allows switching between the three light modes and can also be used during therapy.

Wireless use

The mask is wireless, which significantly increases comfort, allows free movement, and at the same time reduces the risk of connector damage. The wireless solution makes the therapy a comfortable home ritual without limitations.

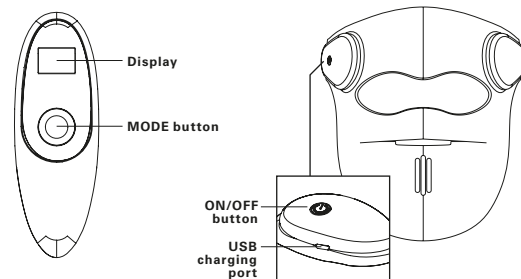2.3. EXPOSURE TIME AND DURATION OF ONE THERAPY SESSION

The LED 2.0 Biohack Mask is set so that one automatic cycle lasts 10 minutes. After this time, the mask switches off and the light stops shining.

Recommended length of one therapy session

To achieve optimal results, we recommend a total duration of one therapy session of 20 minutes. This means that after the first 10-minute cycle ends, you turn the mask on again and complete a second 10-minute cycle.

How to work with modes during a 20-minute therapy session

During one 20-minute therapy session, you may use the same mode during both cycles, or switch between cycles.

2.4. FREQUENCY OF USE

First month: 4x per week.

Second month: 3x per week.

From the third month onward: 2x per week long-term.

3. CONTRAINDICATIONS

The use of the LED 2.0 Biohack Mask is not recommended in the following cases:

- pregnancy or breastfeeding without prior consultation with a doctor
- photosensitivity (increased sensitivity to light)
- use of medications or supplements causing photosensitivity (e.g. certain antibiotics, retinoids, isotretinoin, acne medication, hormone therapy, etc.)
- genetic disorders of the eye or skin
- epilepsy or neurological disorders sensitive to light stimuli
- active oncological disease
- recent surgical or invasive procedures in the facial area without a doctor's approval
- active skin infections, open wounds, or pronounced inflammatory symptoms

In the event of excessive redness, burning, irritation, or any other adverse reaction, stop the therapy immediately and consult a specialist.

4. WARNINGS FOR USE

1. Fully charge the mask before first use. Do not use the mask while charging.

2. Do not expose the device to water or other liquids. Neither the mask nor the remote control is waterproof. Avoid direct contact with water and do not use them in a humid environment (e.g. in the shower or bath).

3. Do not expose the mask to impacts or drops. Mechanical damage may affect the functionality of the device.

4. Do not store the mask in extreme temperatures or in a humid environment. Store it in its original packaging in a dry and cool place, away from direct sunlight.

5. Do not place the device on blankets, electric blankets, or heat sources during operation.

6. Use only the original charging cable supplied by the manufacturer. Do not use a fast-charging adapter, as there is a risk of the battery overheating, shortening of the device's lifespan, and damage to the device.


7. The mask is intended for one person only. For hygiene reasons, it is not suitable to share it with other persons.

8. The device is not intended for children. Keep it out of their reach.

9. Unauthorized assembly, disassembly, or repair of the device may lead to damage and loss of warranty. In the event of a malfunction, contact the seller.

10. Use only 2 AAA batteries. Insert the batteries with the correct polarity according to the marking in the device. Do not use damaged or leaking batteries, do not combine new and used batteries or batteries of different types. If you will not use the device for a longer period, remove the batteries. Batteries are not included in the package.

5. PACKAGE CONTENTS

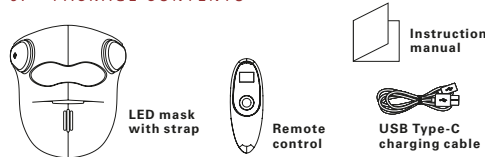6. PRODUCT PARAMETERS

Number of LED diodes: 280 pcs divided into 70 light units
Charging voltage: 5 V
Operating voltage: 3.7 V
Charging current: 890 mA
Input power: 5 W
Output power: 3.7 W
Irradiance: 30–50 mW/cm²
Weight: 320 g
Made of food-grade silicone

Made in Hong Kong exclusively for Envy Therapy
Distributor: Klinika ENVY, s.r.o., Slovakia
envytherapy.com

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