NÁVOD NA POUŽITIE

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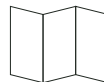
NABÍJECÍ PROUD: 500MA

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TYP BATERIE: 400MAH POLYMER LP

MADE IN HONG KONG EXCLUSIVELY
FOR ENVY THERAPY

DISTRIBUTOR:
KLINIKA ENVY, S.R.O.,
SLOVAKIA

ENVYTHERAPY.COM

CE RoHS

VISIBLY REVITALIZING YOUR SKIN,
THE DEVICE ENHANCES ITS ELASTICITY,
REDUCES PUFFINESS, AND ACCENTUATES
FACIAL CONTOURS.

THE ION TECHNOLOGY ACTIVATES SKIN
CELLS, PROMOTES MICROCIRCULATION.

Envy Therapy Sculpt & Contour ION Stick

RADICAL
SKINCARE SYSTEM

[EN]

Familiarize yourself with the individual functions of the device and the usage warnings provided.

Carefully read and understand all instructions before operating the device.

Ensure the device is fully charged before initial use.


Please note that the device is intended for cosmetic purposes and is not a medical tool.

The effectiveness of the device may vary depending on skin type and the type of cosmetic products utilized.

TABLE OF CONTENTS

1. DEVICE INFORMATION
 - 1.1. MAIN BENEFITS
 - 1.2. MODES DESCRIPTION
2. USAGE INSTRUCTIONS
 - 2.1. HOW TO USE
 - 2.2. DURATION OF USE
 - 2.3. FREQUENCY OF USE
3. WARNINGS
4. PACKAGE CONTENT
5. DEVICE SPECIFICATION

1. DEVICE INFORMATION

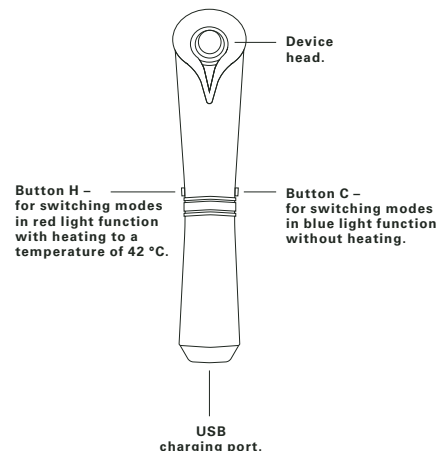1.1. MAIN BENEFITS

Visibly revitalizing your skin, the device enhances its elasticity, reduces puffiness, and accentuates facial contours.

The ION technology activates skin cells, promotes microcirculation, and facilitates the absorption of skincare products due to the presence of positive and negative ions.

Say goodbye to wrinkles and hello to firmer skin as the targeted vibrations awaken your complexion, diminish signs of fatigue, and enhance facial contours.

The unique design ensures seamless skin contact from every angle.


The device is designed for treating the entire face.

Its ergonomic and cordless design offers convenience.

The device can be charged using a USB cable, ensuring safety and environmental friendliness.


1.2. MODES DESCRIPTION

BUTTON H – RED LIGHT
WITH TEMPERATURE OF 42°C

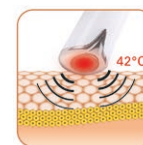
Massage Mode

Upon pressing the H button for the first time, the red light, heating to 42°C, and vibration activate. This massage mode aids in reducing puffiness and enhancing facial contours. Targeted vibrations rejuvenate the skin and diminish signs of fatigue. Use with face cream or oil.

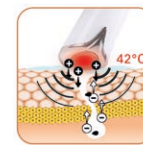Massage Mode + ION Mesotherapy Mode

Pressing the H button for the second time activates the red light, 42°C heating, vibration, and ion mode. The massage experience is further enhanced by the presence of positive and negative ions, promoting better absorption of skincare products into the skin. Suitable for use with a serum or ampoule.

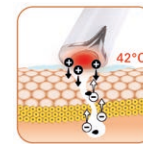
ION Mesotherapy Mode

Upon pressing the H button for the third time, the red light and ion mode engage. Utilize the vibration-free function to enhance the absorption of cosmetic products. Recommended for use with a serum or ampoule.


BUTTON C - BLUE LIGHT WITHOUT HEATING

For individuals with sensitive skin, skin prone to sensitive capillaries, or those preferring a cooler treatment, we recommend utilizing the massage mode without heating.

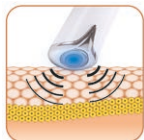
Massage Mode

Press the C button once to activate the blue light and vibration. This massage mode assists in reducing puffiness and improving facial contours. Targeted vibrations revitalize the skin and diminish tiredness. Use alongside face cream or oil.

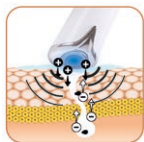Massage Mode + ION Mesotherapy Mode

Pressing the C button for the second time triggers the blue light, vibration, and ion mode. The massage experience is enriched by positive and negative ions, enhancing the absorption of cosmetic products. Use the mode with a serum or ampoule.

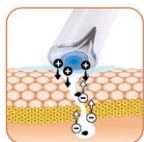
ION Mesotherapy Mode

Upon pressing the C button for the third time, the blue light and ion mode activate. Utilize the vibration-free function to enhance the absorption of cosmetic products. Recommended for use with a serum or ampoule.


2. USAGE INSTRUCTIONS

2.1. HOW TO USE

Apply the product to cleansed skin based on the selected mode (refer to modes description) and follow the illustrated instructions below. For optimal absorption of active ingredients, we recommend applying the product just BEFORE using the device. It is advisable to apply the product in sections.

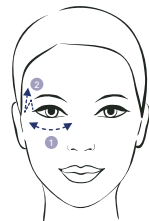Method for applying
to the eye area.

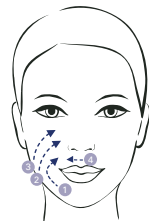Method for applying
to the face.

After use, manually switch off the device if it has not automatically turned off, and wipe the head with a clean cloth. Avoid exposing the device to running water, as it is not waterproof.

2.2. DURATION OF USE

Massage Mode: 5 to 10 minutes.

ION Mesotherapy Mode: 3 to 5 minutes.

The device will automatically shut off after 10 minutes.

2.3. FREQUENCY OF USE

You may use the device daily.

DO NOT USE IF YOU BELONG TO ANY OF THE FOLLOWING GROUPS:

- Pregnant or breastfeeding women.
- Individuals sensitive to light or taking photosensitizing medications.
- Individuals with genetic eye or skin conditions (e.g., lupus erythematosus, photosensitive eczema, albinism, etc.).
- Individuals prone to epileptic seizures.
- Individuals with pacemakers.

3. WARNINGS

1. Ensure the device is fully charged before the first use and avoid using while charging.
2. Refrain from cleaning, washing, or using the device during charging. Keep away from direct water contact.
3. The device features a heating function with an approximate 42°C temperature. Sensitivity to temperature varies among individuals.
4. Discontinue use and seek advice from a beautician or doctor if excessive redness or adverse reactions occur.
5. If you suspect any issues with the device, discontinue use and contact the seller.
6. The product is not waterproof; avoid contact with water.
7. Handle the device with care; do not drop or bump it.
8. Store the device in its original packaging in a cool, dry place away from extreme temperatures and moisture. Keep it out of direct sunlight.
9. Use only the original charging cable.
10. Be careful when using and storing the device. Keep it out of the reach of children.
11. The device is not intended for consumption.
12. Do not bring it into contact with liquids, metal products, or corrosive chemicals.
13. Unauthorized assembly, disassembly, repair, and modification of the device can lead to malfunctions that cannot be repaired or claimed. Do not attempt to repair it yourself and always contact the distributor.

Following these instructions will help prevent damage to the device and ensure your safety while using it.


4. PACKAGE CONTENT

Sculpt & Contour
ION Stick

USB
charging cable

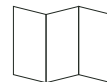instruction
manual

5. DEVICE SPECIFICATION

MODEL OF THE DEVICE: SY-026

WEIGHT: 60 G

CHARGING VOLTAGE: 5V

CHARGING CURRENT: 500MA

RATED POWER: 3W

BATTERY TYPE: 400MAH POLYMER LP

MADE IN HONG KONG EXCLUSIVELY
FOR ENVY THERAPY

DISTRIBUTOR:
KLINIKA ENVY, S.R.O.,
SLOVAKIA

ENVYTHERAPY.COM

CE RoHS