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THE ULTRA-INTELLIGENT HI-TECH DEVICE FOR ALL SKIN TYPES WORKS ON THE PRINCIPLE OF ULTRASONIC WAVES IN 3 DIFFERENT MODES:

- 1—ULTRASONIC CLEANSING MODE
- 2—THE GALVANIC MODE
- 3—EMS STIMULATION OF FACIAL MUSCLES

Envy Therapy
Ultrasonic
Face Therapy
3-in-1

RADICAL
SKINCARE SYSTEM

This user manual explains the individual functions of the device as well as warnings for safe use.

Before using the product, please read and follow the user manual in detail.

Charge the device before using it for the first time.

The device is not a medical, but a cosmetic tool.

The skin type and the type of cosmetic products used affect the effectiveness of the device.

The product has been tested by the Korean Institute of Dermatological Sciences.

CONTENTS




- 1. INFORMATION ABOUT THE TOOL
 - 1.1. PRODUCT SPECIFICATION
 - 1.2. PARTS OF THE DEVICE
 - 1.3. NAMES OF INDIVIDUAL PARTS
- 2. USE
 - 2.1. CHARGING
 - 2.2. MANUAL
- 3. WARNINGS

1. INFORMATION ABOUT THE TOOL

1.1. PRODUCT SPECIFICATION


NAME: ULTRASONIC FACE THERAPY 3-IN-1
 BATTERY: 500 mAh@3.7 V
 SIZE: 24x63x180 mm (WXLXH)
 WEIGHT: 120 g
 MATERIAL: ABS + STAINLESS STEEL
 PRODUCTION YEAR: 2024
 COUNTRY OF ORIGIN: PRODUCED IN HONG KONG
 EXCLUSIVELY FOR ENVY THERAPY
 MAIN FUNCTIONS: ULTRASONIC DEEP CLEANSING,
 GALVANIC PENETRATION OF ACTIVES, EMS

1.2. PARTS OF THE DEVICE

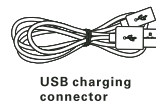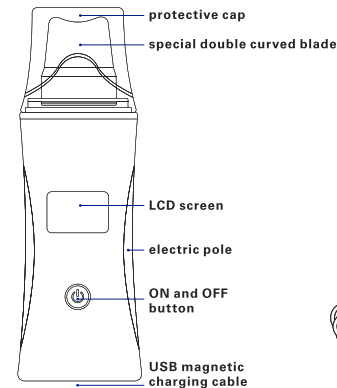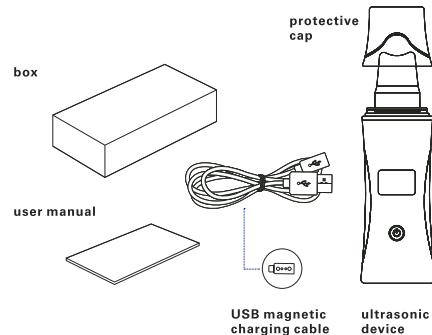2. USE

2.1. CHARGING

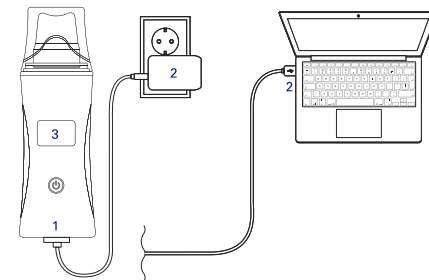Charge the tool before using it for the first time. The condition of the battery and charging time may vary. Maximum charging time is 2 hours.

- 1. Connect the microUSB magnetic cable to the connector on the bottom of the tool.
- 2. Connect the other part of the cable to a 5V 1A plug or a computer.
- 3. The battery symbol on the LCD screen flashes during charging. Do not use the product while charging.

WARNINGS

If the battery symbol flashes on the LCD screen during use, it means the battery is dying and the tool needs to be charged.

Do not use the product while charging.

The battery symbol on the LCD screen flashes during charging.

When the battery symbol stops flashing, the tool is fully charged and ready for use.

Once fully charged, the tool can be used for up to one hour. The time may vary depending on the operating mode used.

2.2. MANUAL

Before using the product, please read and follow this manual in detail.

The tool has an automatic shut-off function after 5 minutes of operation. If you need to use the device for more than 5 minutes, please press the power button again and select the mode you need. When using an EMS mode, we recommend starting with EMS/L and then changing to EMS/H. Alternatively, leave on EMS/L depending on your subjective feelings during use. Always clean the instrument with a disinfectant after use.

A) ULTRASONIC DEEP CLEANSING (SONIC CLEANSING)

High frequency vibrations help cleanse the pores from dead skin cells and bacteria.

Before using the tool, remove make-up, clean your skin and moisturize it. Make sure that the treated area is always wet. You can use plain water or a mist. Use a tampon or spray to moisten the treated area continuously. However, do not use products that contain alcohol.

Turn on the device and select **Sonic Cleansing**. Gently move the blade vertically across the skin (Fig. 1).

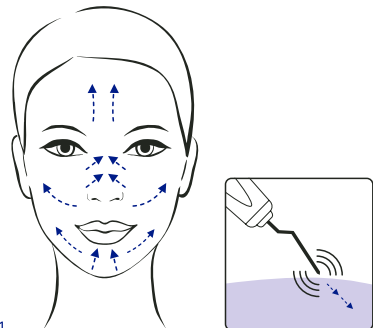
FIG. 1

B) ULTRASONIC DEEP CLEANSING ION+ (SONIC CLEANSING ION+)

High frequency vibrations help cleanse the pores from dead skin cells and bacteria. In addition, positively charged ions help remove cosmetic residual and synthetic impurities (dust, smog, dirt...) from the skin.

Before using the tool, remove make-up, clean your skin and moisturize it. Make sure that the treated area is always wet. You can use plain water or a mist. Use a tampon or spray to moisten the treated area continuously. However, do not use products that contain alcohol.

Turn on the device and select **Sonic Cleansing ION+**. Gently move the blade vertically across the skin (Fig. 1).

FREQUENCY AND PROTOCOL FOR USING OF ULTRASONIC DEEP CLEANSING MODE

- OILY SKIN – TWICE A WEEK:
5 MINUTES SONIC CLEANSING
5 MINUTES SONIC CLEANSING ION+
- NORMAL SKIN – ONCE A WEEK:
5 MINUTES SONIC CLEANSING
5 MINUTES SONIC CLEANSING ION+
- DRY AND SENSITIVE SKIN – ONCE A WEEK:
5 MINUTES SONIC CLEANSING

C) GALVANIC MODE ION-

Galvanic ion booster helps the actives from skincare products penetrate deeper into the skin.

Before using the instrument, remove make-up and clean your skin. Apply serum, concentrate or ampoule according to the needs of your skin. For maximum penetration effect of active substances, we recommend applying serums, concentrates and ampoules just BEFORE using galvanic mode. Ideally, apply the product and use galvanic mode in sections. For example, apply serum to the forehead, use a galvanic mode. Continue similarly in the sections: right cheek, left cheek, chin, neck, décolleté.

Select the galvanic mode (**Galvanic ION-**) on the tool and hold it by the electric poles on the sides. Then, with the back of the blade (not its tip), gently press the skincare product into the skin, ideally in sections (Fig. 2).

FREQUENCY AND PROTOCOL FOR USING OF THE GALVANIC MODE

- 3-5 MINUTES 2-3 TIMES A WEEK AS NEEDED

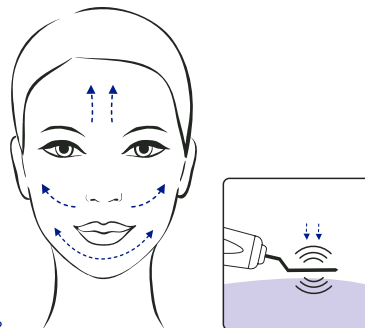
FIG. 2

D) ELECTRIC FACE MUSCLE STIMULATION (EMS/LOW)

The EMS (electric muscle stimulation) lifts your face contours and improves skin tightness.

Before using the instrument, remove make-up and clean your skin. Apply the serum, concentrate, ampoule or cream according to your skin needs. Make sure to use a water-based, not an oil-based product. Wait until the product soaks into the skin. When using EMS mode onto wet skin, an increased tingling sensation may occur.

Select the **EMS/L** mode and hold the tool by the electric poles on the sides. Then, with the back of the blade (not its tip), gently massage the skin (Fig. 2). It is normal to feel a tingling sensation while using this mode.

E) ELECTRIC FACE MUSCLE STIMULATION (EMS/HIGH)

Higher frequency of EMS (electric muscle stimulation) for even better result in face contours lifting and skin tightness.

Before using the tool, remove make-up and clean your skin. Apply the serum, concentrate, ampoule or cream according to your skin needs. Make sure to use a water-based, not an oil-based product. Wait until the product soaks into the skin. When using EMS mode onto wet skin, an increased tingling sensation may occur.

Select the **EMS/H** mode and hold the tool by the electric poles on the sides. Then, with the back of the blade (not its tip), gently massage the skin (Fig. 2). It is normal to feel a tingling sensation while using this mode.

If the frequency of this mode is too high and causes you discomfort, switch to EMS/L, a lower frequency.

FREQUENCY AND PROTOCOL FOR USING THE EMS MODE

- INTENSIVE EMS TREATMENT:
5 MINUTES 4 TIMES A WEEK FOR ONE MONTH
- MAINTENANCE EMS TREATMENT:
5 MINUTES 3 TIMES A WEEK FOR ONE MONTH
- PREVENTIVE EMS TREATMENT:
5 MINUTES 2-3 TIMES A WEEK AS NEEDED

3. WARNINGS

Before using and storing the device, it is necessary to read the following warnings.

The water resistance level of the product is IP X-5. Make sure no water or other liquid enters the device.

Ultrasound treatment is not recommended for persons with a pacemaker or for persons suffering from epilepsy or severe heart disease. Do not use on skin with ongoing acute inflammation or fresh permanent make-up.

Due to the complexity of the issue and serious ethical issues, clinical trials of the device have not been performed in pregnant women and we therefore do not recommend using it during this period.

If excessive redness or other adverse reactions occur, stop using the device immediately and consult a beautician or doctor.

If you suspect that something is wrong with the device, stop using it immediately and contact the seller.

Do not throw or thump the product.

Do not store at temperatures too high or too low or in humid environments. Store in the original box in a dry and cool place. Do not expose to direct sunlight.

Do not place the tool on regular or electric blankets.

Use only the original charging cable that came with this product.

Be careful when using and storing the device. The device is not suitable for use by children, so keep it out of their reach.

The device is not intended for consumption.

Do not expose the product or the charging cable to liquids or other metal products.

Arbitrary assembly, disassembly, repair and modification of the device may lead to malfunction that will not be possible to repair or it will not be possible to process a complaint. Do not attempt to repair the device yourself and always contact the seller.